# F\_Arco



**Fedesa** 

# **F\_Índice**



# **Fedesa**

Since 1970, two family generations and a large professional team have devoted all their passion and knowledge towards creating the best experience in dental care for dentists and their patients.

We design and manufacture high performance dental units in Spain, in which every component has been developed to be comfortable, efficient and long-lasting.



**3\_**Finishes

**1\_**About us



Your essential teammate, always by your side

FEDESA is a dentist's essential teammate - equipment you can rely on for doing your job and a team that shares the same passion for dental health.

A dental chair always goes alongside a dentist. It's their best ally for inspiring trust and peace of mind in patients.

At FEDESA, we are part of your team, and we are always at your side so that nothing gets in the way of your patients' treatment.

# Fedesa High performance, maximum reliability

FEDESA has been the market leader in high performance dental units in Spain for over 50 years.

With our experience in design, processes and manufacturing technology, we are able to make high performance equipment accessible to all kinds of clinics.

All our work is focussed on offering yu the best experience, service and durability.

Dentist take centre stage and are the core of all that we do. They have a very evident vocation that defines them as individuals.

This vocation is shared by the whole FEDESA team. A family united by a passion to help professionals in the dental sector.



# A great, leading company that casts its legacy towards the future

FEDESA has been the undisputed leader of the Spanish market for decades. As well as leading in terms of market share, we are a great company because of our history and the strength of our human resources.

We do not settle for what has already been done. Rather, we looktowards the future with an ambition for continual improvement, innovation and sustainability.



# High performance dental units for the most discerning professionals

The most discerning dentists choose FEDESA because our high performance dental units are designed and manufactured to be the essential teammate that makes your job easier for many years to come.

We guarantee the performance and durability of our units as a result of our rigorous production processes, quality control and external certifications.



# A committed, passionate team that is always there for dentists

The FEDESA team shares the same values that it has transmitted from generation to generation.


All our work aims to make a dentist's job easier.





# For the most discerning professionals

We design and manufacture high performance dental units in Spain, in which every component has been developed to be comfortable, efficient and long-lasting.

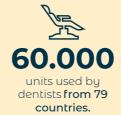





anywhere in Europe



production and development centre in **Madrid** measuring

Dentists treat patients with serious oral health problems. They need to diagnose, treat and follow up with patients until they have recovered.

Dentists also conduct preventative treatments to ensure their patients' oral health.

Just like their patients, dentists need quick responses and the ability to react in adverse circumstances.

We conduct maintenance and follow-up to ensure maximum performance for your equipment.

FEDESA is a dentist's essential teammate, always at their side, sharing the same passion for oral health.



# ¿What is the Arco model?

a

10

The Arco dental unit, Fedesa's latest product, appeared on the market in 2012. The unit, which fits in to the medium-to-high range market niche, brings together all the standard resources needed by modern dentists.

All the components of the unit - dental chair, cuspidor unit, main delivery system and auxiliary tray - boast the **most highly advanced systems** and technology.



hanging hoses and with a mobile unit.



Dental Unit

Fedesa

# F\_Arco





#### **Arco Premium**

The Arco Premium unit incorporates in its basic configuration the following equipment: an independent control panel for the dental unit and for the control of the micromotor and scaler, LED light on all the instruments, endo induction micromotor, ultrasonic scaler with LED light and a LED type ilumination lamp.

Page 14-15

#### **Arco Lux**

The most standard configuration, the unit's safety systems ensure incident-free use. Safety system applied to the chair lowering function, to the backrest and to the assistant's workstation.

Incorporates multi-purpose hoses with light, an electric brushless or induction micromotor with light and ultrasonic scaler.

Page 16-19



### Arco Eco

The **Eco version** is a pneumatic unit with the cuspidor and tray attached to the structure of the chair.

The assistant's position features a simple tray with two holders and the main tray includes pneumatic instruments with an ultrasonic scaler.

Page 20-23



### **Arco Orthodontics**

The Arco dental unit also offers a special configuration for Orthodontic treatments. All the instruments required for these purposes can be found in the assistant's work position.

This version also features an instrument tray that is fitted on a double-articulated arm that provides a wide range of movements.


Page 24-25

# F\_Arco

### **Arco Premium**

This model incorporates in its basic configuration the following equipment: an independent control panel for the dental unit and for the control of the micromotor and scaler, LED light on all the instruments, endo induction micromotor, ultrasonic scaler with LED light and a LED type ilumination lamp.



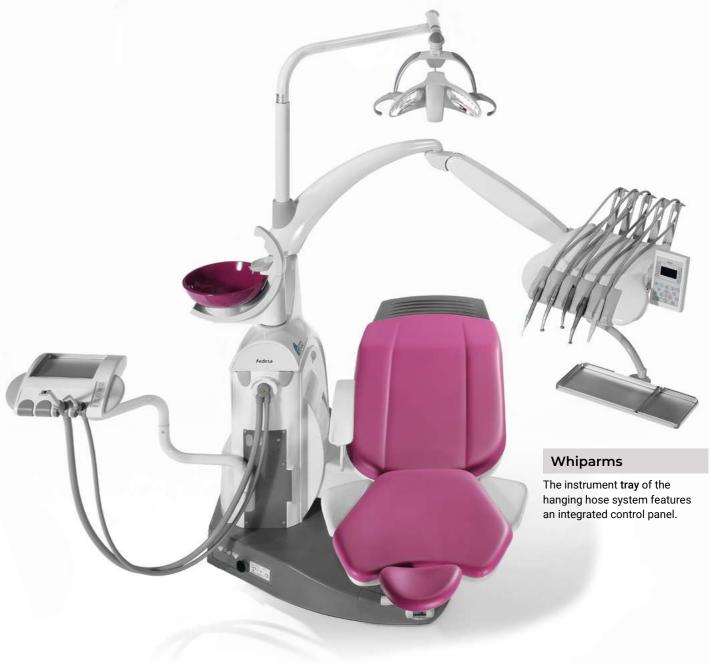



#### Headrest

Height adjustable and double articulation **headrest**. Its ergonomics allow patients to place their heads in the ideal position for the doctor to work.

14



## **Fedesa**

### **Arco Lux**

The most standard configuration, the unit's safety systems ensure incident-free use. Safety system applied to the chair lowering function, to the backrest and to the assistant's workstation. Incorporates multi-purpose hoses with light, an electric brushless or induction micromotor with light and ultrasonic scaler.



#### Foot pedal

The unit features a foot pedal to control the speed of the instruments. Optionally, foot pedal can also control movement of the chair.

#### Cuspidor of unit

Is designed to include additional elements on assistant tray such as polymerization lamp, pneumatic hose, syringe, etc...


Optionally can be installed all kind of aspiration systems (wet, dry, amalgam separator, centralized system, etc...).

Cuspidor is also ready to get decontamination system in handpieces, pressurized water tank with selector, hot water, automatic cleaning of aspiration hoses, etc...

**C**rco

# F\_Arco



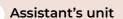

Includes a tray with four holders on a double-articulated arm that allows the assistant to place the unit in a wide range of positions.



#### Water bowl

The water bowl can spin 180° for the benefit of the assistant, making it easier for patients to rinse their mouths. The bowl and faucet of the water supply system enjoy an excellent level of mobility (with a working radius

of + - 90°) to adapt to the various needs of patients.

19 18

## **Fedesa**

# F\_Arco

### Arco Eco

The **Eco version** is a pneumatic unit with the cuspidor and tray attached to the structure of the chair. The assistant's position features a simple tray with two holders (four can be provided on request) and the main tray includes pneumatic instruments with an ultrasonic scaler.



(drco


#### Lockable Whiparms

To avoid the stress when handling the instruments. Reduced height of the arm to avoid contact with other nearby elements. Long hoses.



#### The backrest

Optionally, the entire Arco dental unit range can be manufactured with two types of backrest: the standard model that is supplied with the unit and a model with a more pronounced groove.

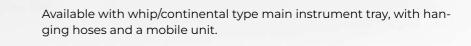




# **Fedesa**



24 25

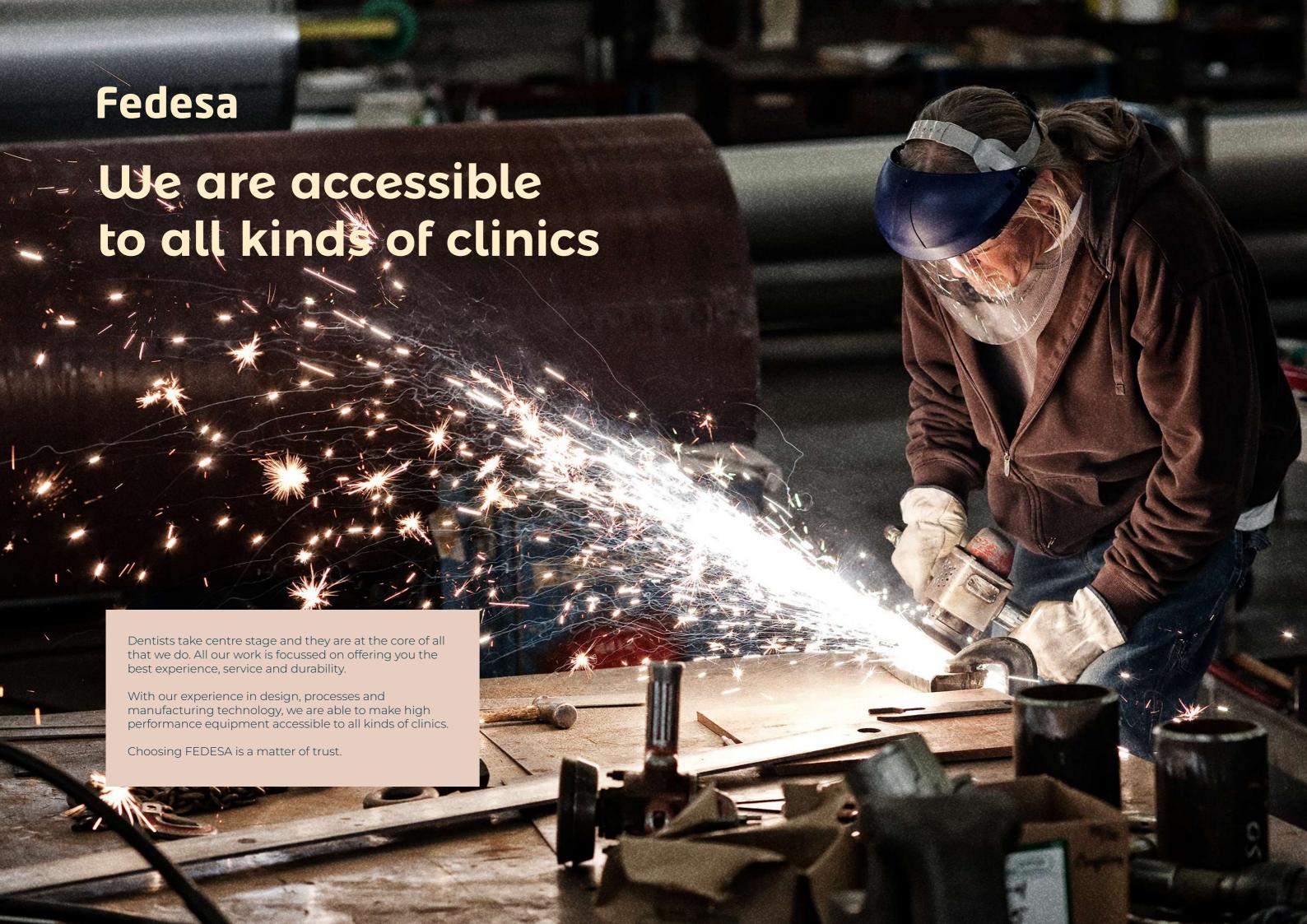

# Finishes

# F\_Arco

#### Textile colour finishes that combine with the dental unit

We have selected a wide range of colours for standard use to integrate with modern facilities.



#### Upholstery

FEDESA provides upholstery on all its dental units that add to the design of its equipment, improving their convenience and comfort.


Low maintenance flame retardant upholstery without seams, washable and excellent resistance to use.


The extensive catalogue of upholstery colours also comes with viscoelastic upholstery (memory foam).

#### **Measures of the dental unit**

The measurements indicated on the drawing will help you place your equipment in your surgery.

#### Dimensions in mm.

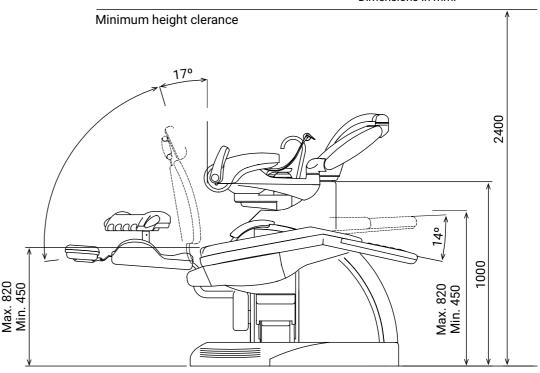

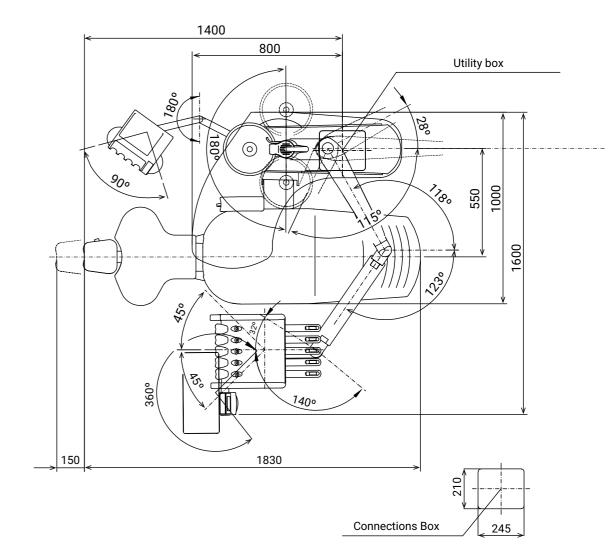

# Options

# F Arco

#### **Ilumination lamp**

All FEDESA dental units include lighting equipment that can be changed as an option.

EDI Lamp (25,000 lux halogen lamp). MAIA Lamp (35,000 lux LED lamp). ALYA Lamp (50,000 lux LED lamp).







Maia Model



Alya Model

#### **Control pedals**

There are two types: pneumatic pedals with instrument rpm control, water and chip blower selector, and a multifunction pedal for instrument control and to adjust the chair.



Progresive Pneumatic Model



Multifunction Model



Pneumatic Model

#### **Footrest**

All seats on FEDESA units can be supplied with footrests, either with transparent plastic protections or independent polyurethane footrests.





#### **Monitor**

Slim, lightweight, with a 17" screen and adaptable to the various positions required for use by the doctor and patient.



#### Intraoral camera

High-resolution camera that can be added to the main tray and to the auxiliary tray. Small and very light, touch activation and detachable.



### Mobile Units

There are various ranges of mobile units: **Premium**, **Lux** and **Eco**. The Premium range features an induction micromotor integrated into the control panel. The Lux range has an electric brush micromotor and the Eco range features pneumatic instrumentation.



Arco Premium Model



Arco Lux Model



Arco Eco Model

### Stools

FEDESA has three types of stools (the **Astral**, **Pony** and **Sory** models) available with or without backrests and with a footrest ring. Surround lumbar support.

They are manufactured in stainless steel and feature a five-wheel base, as well as gas piston height adjustment.

A wide range of colours and textures in the upholstery catalogue.



Astral - A Model

Arco Model




Pony Model with backrest



Pony Model without backrest



Sory Model

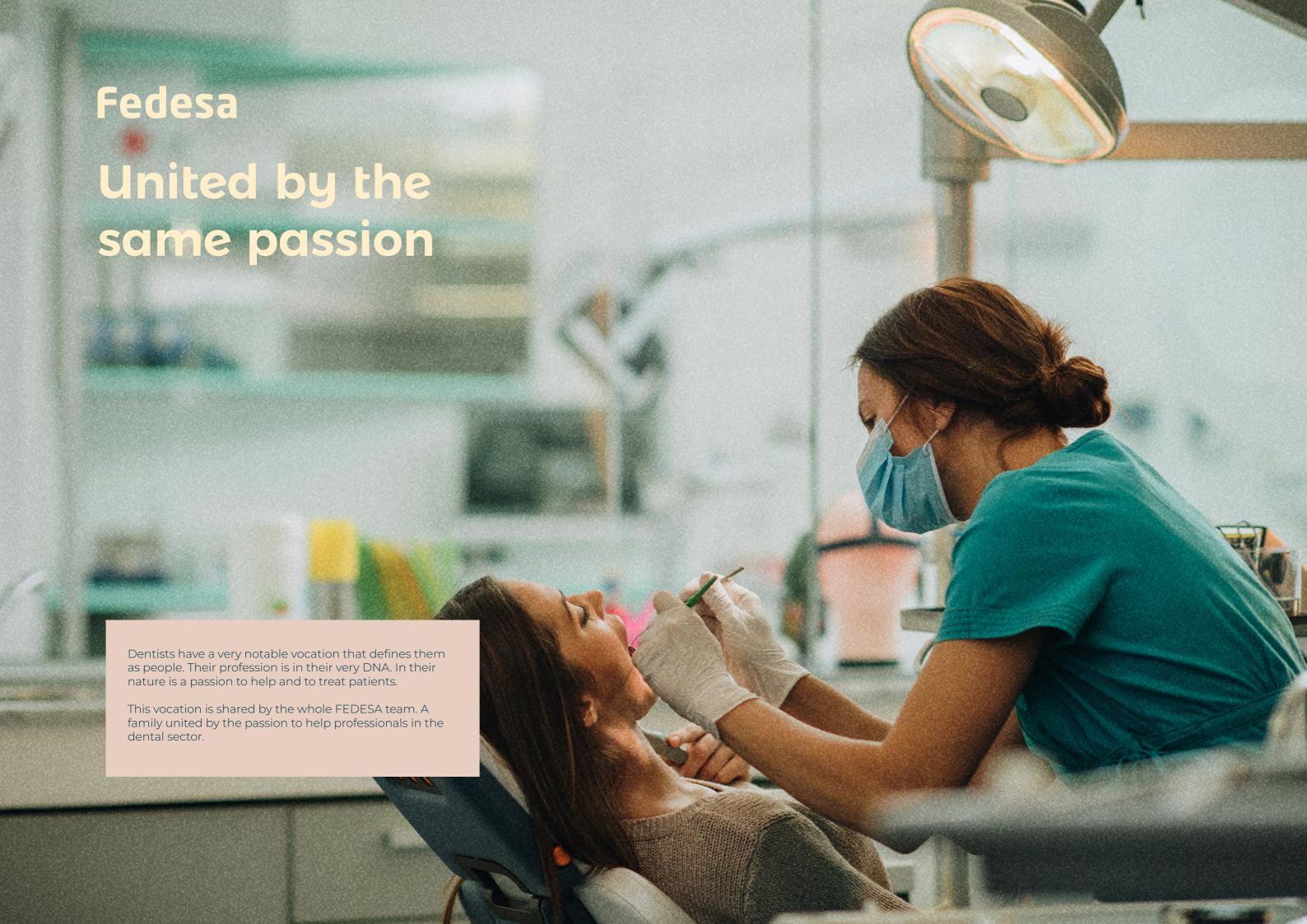



#### **Manufacturing dental units since 1970**

From the origins of **FEDESA** back in the 1970s and right up to the present day, we have dedicated our efforts to meet the needs of dentists.


We are experts in designing, manufacturing, and marketing dental units. Our extensive catalogue of solutions caters for all the needs of today's dentists.

Our products meet all national and international regulations and certifications that ensure the quality of all our products. A wide network of agents approved by our brand for the domestic market, as in the rest of the world, collaborate with us to offer the best customer sevice.



FEDESA is a global company with an unbeatable history in quality, technology and service.

It has an extensive network of distributors all over the world and an excellent professional team that is committed to the company.

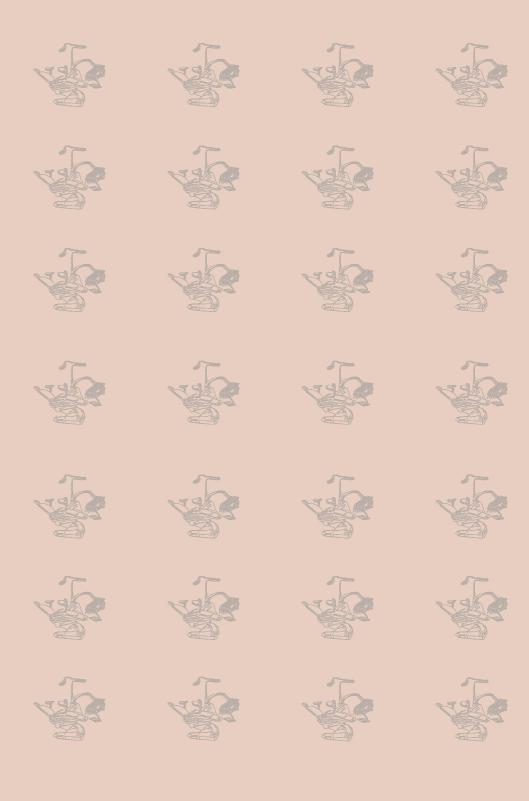

# **Fedesa**

#### NUEVA FEDESA, S.A.

Avda. Madrid, n° 45. 28500 Arganda del Rey (Madrid) España Tel.: 34 - 91 871 23 83. Fax: 34 - 91 871 64 88 fedesa@fedesa.com - www.fedesa.com

FEDESA se reserva el derecho de introducir cualquier cambio en beneficio del producto CT9A23-II. (1000 EJ. 20-03-019) Diseño: De Ledesma. Imprime: Gráficas Aries. Depósito legal: M-7XXXXX

